



Software Updates

Release of Noxturnal 7.0

At Nox Medical, we are committed to continuous innovation and excellence. With the release of our latest Noxturnal 7.0 software update, you will discover performance improvements and new advanced features, all crafted to enhance your experience with Nox Medical's sleep diagnostic technology.

Main updates include:

- AASM compliance for reporting of 3 and 4%
- User authentication
- Audit logs in Noxturnal
- Important stability improvements and enhancements

Download Noxturnal 7.0

<https://noxmedical.com/downloads>

Download Release Bulletin

<https://noxmedical.com/release-bulletin-hub/>

Recommended Practices

Please see the document provided alongside

Software Download

To download the Noxturnal 7.0 software update, please visit our website.

Release Bulletin

Within our Release Bulletin Hub, you can access a comprehensive document detailing the advanced features and enhancements included in the Noxturnal 7.0 software update.

Upgrading to Noxturnal 7.0 - Recommended Practices

We have compiled a document outlining best practices for identifying which settings will be maintained during an upgrade and which will reset to defaults when upgrading to Noxturnal 7.0.


Noxturnal Upgrade

- Recommended practises

Which settings will be kept during upgrade and which settings will be restored to factory default


	Admin Mode (default)		User Mode (custom)	
	Not affected	Reset to default	Not affected	Reset to default
Analysis Protocols	If saved under new name	YES	✓	
Device Profiles	If saved under new name	YES	✓	
Event Types		YES	✓	
Event Groups		YES	✓	
Keyboard shortcuts		YES	✓	
Colours		YES	✓	
Workspace Layouts	If saved under new name		✓	
Battery Profile	If saved under new name	YES	✓	

The difference between user database and default database

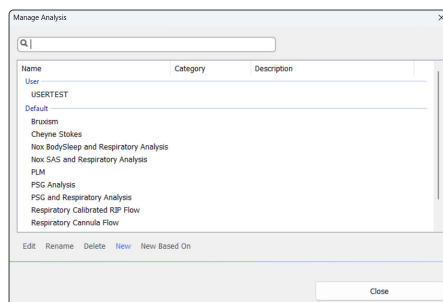


Analysis Protocols

Analysis listed in the user section will be maintained.
Analysis in the default section will be reset to default

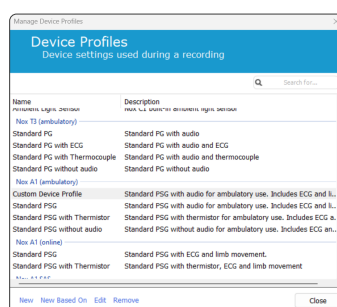

Device Profiles

Default device profiles will be reset to default during upgrades.
Custom device profiles will be maintained and can be detected as they can be edited and removed

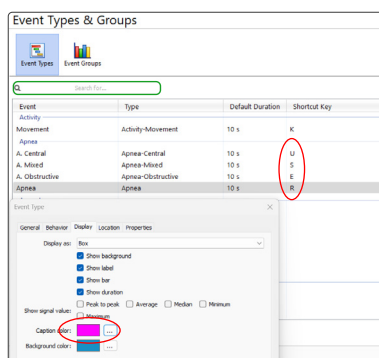

Event Types and Groups

Shortcuts and colors created in user mode will be maintained during upgrade. If created or changed in admin mode, shortcuts and colors will be reset to default.

When updating the Noxturnal language pack, there is no requirement to uninstall the previous version. The new version will automatically overwrite the existing one